## Benchmark



- Temperature control
  from 0° to 100°C

  Thocks.
- Precision contoured wells, for uniform thermal transfer
- Two models: heat/cool or heat only

## **Technical Data:**

Temp. Range (HC): Temp. Range (H): Heating Ramp Rate: Cooling Ramp Rate: Temp. Accuracy: Speed Range: Motion / Orbit: Capacity: Dimensions:

Ambient +5° to 100°C 5.0°C / min. 2.5°C / min. (HC only) 200 to 1500 rpm, or "off" Horizontally circular / 2mm Per selected block 8.25 x 11.5 x 8.25 in. 21 x 29 x 21 cm 18lb. / 8kg 120V, 60Hz or 230V, 50-60Hz

6 x 50ml

12 x 15ml

Ambient -20° to 100°C

Weight: Electrical:

useful for a variety of molecular biology applications including: - Denaturation of DNA, RNA and proteins - Reverse transcription - Lipid extractions

The MULTITHERM shaker is a temperature controlled vortexer,

- Yeast and bacteria cultivation


- Restriction digests

- Gentle thawing of biological materials

- DNA/plasmid isolation

 $\epsilon$ 

Despite its modest 8 x 11 in. footprint, all popular tubes and plates are compatible, including as many as six 50ml tubes. Speed, time and temperature settings are continuously visible on the LCD, simultaneously showing both actual and selected values. Integral over-temperature control ensures long life, safety and sample integrity.

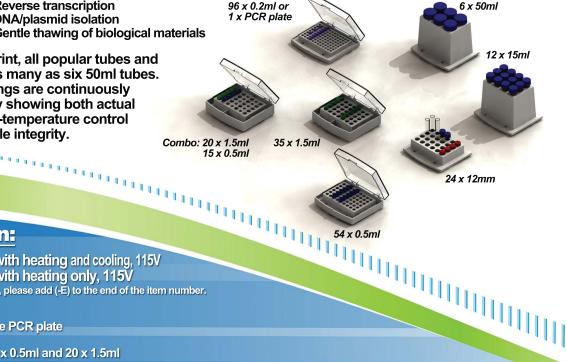

## **Ordering Information:**

H5000-HC MultiTherm Shaker with heating and cooling, 115V H5000-H MultiTherm Shaker with heating only, 115V

To order in 230V, please add (-E) to the end of the item number.

Please order blocks separately:

H5000-02 Block, 96 x 0.2ml or one PCR plate

H5000-05 Block, 54 x 0.5ml

H5000-CMB Block, combination 15 x 0.5ml and 20 x 1.5ml

Block, 35 x 1.5ml Block, 35 x 2.0ml H5000-15 H5000-20

H5000-12\* Block, 24 x 12mm (200 to 1200rpm / Ambient -10° to 100°C) H5000-150\* Block, 12 x 15ml (200 to 750rpm / Ambient -10° to 100°C) Block, 6 x 50ml (200 to 750rpm / Ambient -10° to 100°C) H5000-500\*

\* Note: The larger mass of this block results in a reduced speed and temperature range.

**Benchmark Scientific Inc.** PO Box 709 - Edison, NJ 08817 Phone: 908-769-5555

Fax: 908-222-1864

Web: www.BenchmarkScientific.com Email: Info@BenchmarkScientific.com