





# Compact Low-Speed Centrifuge with a Modern Touch-Wheel Interface Ideal for General Lab Applications

The FC5707 is a reliable compact low-speed general purpose centrifuge which includes a pre-installed  $8 \times 15$  mL rotor. The rotor can process both round and conical bottom tubes (including common blood tubes), and can also be adapted easily with accessories to process smaller tubes. Designed with convenient operation in mind, the FC5707 features a uniquely designed splash-proof front panel with Touch-Wheel control on an intuitive interface.

#### **Standard Features:**

- Modern Touch-Wheel Control on a Sleek Splash-Proof Interface—Designed for ease of operation, the FC5707 features convenient control on a protective splash-proof interface, which allows for responsive operation—even when wearing gloves.
- The Compact Design Saves Valuable Benchtop Space—Enclosing all its functionalities within a small footprint, Frontier<sup>™</sup> Multi Centrifuges save valuable space on the lab bench.
- Easy and Flexible Operation for a Range of Applications—The 8 × 15 mL angle rotor can process both round and conical bottom tubes, including common blood tubes. It can also be adapted for use with smaller tubes such as 7, 5 and 1.5/2 mL when used with accessories.

# Frontier<sup>™</sup> 5707 Multi Centrifuge

Model	FC5707		
Speed Range	200 rpm – 6800 rpm; 10 rpm/set		
Maximum RCF	4445 × g; 10 × g/set		
Maximum Capacity (Rotor)	8 × 15 mL RB* / 4 × 15 mL FA*		
Temperature Range	Air cool		
Running Time	Quick Spin, 10 sec to 99 hr 59 min 59 sec, or continuous		
Noise Level	≤ 60 dB(A)		
Acc/Dec Levels	10		
Allowable Density at Maximum Speed	1.2 g/mL		
Allowable Kinetic Energy	845 Nm		
AC Power	230 V ~ 50/60 Hz		
Voltage Fluctuation	± 10%		
Current Consumption	0.8 – 0.4 A		
Power Consumption	50 W		
Dimensions (W $\times$ D $\times$ H)	277 × 356 × 236 mm		
Net Weight (with rotor)	10 kg		
Shipping Dimensions ( $W \times D \times H$ )	370 × 470 × 340 mm		
Shipping Weight (with rotor)	11.5 kg		
Environment	For indoor use only		

\*RB-round bottom, FA-conical bottom

## **Other Standard Features and Equipment**

Ten acceleration and ten deceleration rates, pre-selection of running time from 10s to 99h 59min or continuous, audible signal at the end of each run with melody options, quick spin for short runs


#### Compliance

- Product Safety: IEC/EN 61010-1; IEC/EN 61010-020
- Electromagnetic Compatibility: IEC/EN 61326-1 Class B, Basic environments
- Compliance Mark: CE

## **Rotor and Accessories**

Order No.		30472305			
	Description	Rotor, Angle, 8 × 15 mL			
	Max. Speed	6800 rpm			
Max. RCF		4445 × g			
Accessories					
Order No.	Description	Vol./Tube	Diameter	Tube type	
Without Adapter	Rotor, Angle, 8 × 15 mL Round bottom tubes/ 4 × 15 mL Conical bottom tubes	15 mL	17 mm	Round bottom/ Conical bottom	
30130889	Adapter, 1 × 7 mL D13.5 mm RB, 2/pk	7 mL	13.5 mm	Round bottom	
30130890	Adapter, 1 × 5 mL D13.5 mm RB, 2/pk	5 mL	13.5 mm	Round bottom	
30130886	Adapter, 1 × 1.5/2.0 mL D11 mm, 6/pk	1.5/2 mL	11 mm	Microtube	

### Dimensions

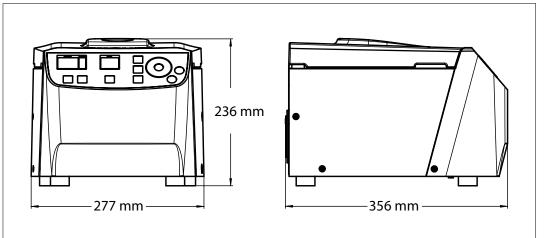

#### **OHAUS Europe GmbH**

lm Langacher 44 8606 Greifensee Switzerland

e-mail: ssc@ohaus.com Tel: 0041 22 567 53 19 e-mail: tsc@ohaus.com Tel: 0041 22 567 53 20

#### www.ohaus.com

The management system governing the manufacture of this product is ISO 9001:2015 certified.



80775550 20181024 © Copyright OHAUS Corporation

