

**AF7015** (art. 904569)

**ABSOLUTE HEPA FILTER KIT**

- 2 absolute HEPA filters
- class EN 1822-5 H14 99.995% for drying air



**Machines**

GW6290DS                      GW7015                      GW7015T

**AS4060** (art. 902897)

**AQUA STOP WATER LEAK KIT**

- prevents accidental water leakage
- detects water leaking from machine and stops it



**Machines**

GW0260-00-00J	GW0260-S0-00J	GW1260-00-000	GW1260-S0-000
GW1260-SC-000	GW4260-00-000	GW4260-S0-000	GW4260-S0-060
GW4260-SC-000	GW4290-00-000	GW4290-S0-000	

**AS4190** (art. 903079)

**AQUA STOP WATER LEAK KIT**

- prevents accidental water leakage
- detects water leaking from machine and stops it



**Machines**

GW4290-00-000                      GW4290-S0-000

## **DRAINP6010** (art. 903465)

---

### **DRAIN PUMP FOR WALL DRAIN**

- additional pump for wall drain
- to be used when the standard floor drain is not available
- hose diameter: 21 mm



### **Machines**

GW6010L	GW6010LX	GW6010ML	GW6010MLX	GW6010T
GW6010TM	GW7010	GW7010M	GW7010T	GW7010TM

## **DRAINP7015** (art. 904265)

---

### **DRAIN PUMP FOR WALL DRAIN**

- additional pump for wall drain
- to be used when the standard floor drain is not available
- hose diameter: 21 mm



### **Machines**

GW7015	GW7015T
--------	---------

**ESVC10** (art. 904961)

**EXHAUST AIR BREAKER 6010-7010**

- Chimney vent disconnector, to be applied to the smoke vent in the upper part of the machines
- It avoids the transmission of pressure variations to the machine and the consequent alarm signaling
- The special shape is designed to prevent condensation from returning inside the machine
- Suitable for 6010/7010 series machines



**Machines**

GW6010L	GW6010LX	GW6010ML	GW6010MLX	GW6010T
GW6010TM	GW7010	GW7010M	GW7010T	GW7010TM

**ESVC15** (art. 904962)

**EXHAUST AIR BREAKER 7015**

- Chimney vent disconnector, to be applied to the smoke vent in the upper part of the machines
- It avoids the transmission of pressure variations to the machine and the consequent alarm signaling
- The special shape is designed to prevent condensation from returning inside the machine
- Suitable for 7015 series machines



**Machines**

GW7015	GW7015T
--------	---------

**ESVC6290** (art. 904963)

**EXHAUST AIR BREAKER 6290**

- Chimney vent disconnectors, to be applied to the smoke vent in the upper part of the machines
- They avoid the transmission of pressure variations to the machine and the consequent alarm signaling
- The special shape is designed to prevent condensation from returning inside the machine
- Suitable for 6290 series machines



**Machines**

GW6290DS                      GW6290DSG

**IC5090** (art. 901737)

**QUALITY WASHING CHECK CONDUCTIVITY SENSOR**

- allows to detect impurities in the washing lye
- activates further rinse phases for achieving the perfect rinse



**Machines**

GW3060BX	GW3060BXA	GW3060BXC	GW3060BXCA	GW3060P
GW3060PA	GW3060S	GW3060SA	GW3060SC	GW3060SCA
GW4090	GW4090A	GW4090BXC	GW4090BXCA	GW4090BXP
GW4090C	GW4090CA	GW4090P	GW4090PA	

**IC6010** (art. 903527)

**QUALITY WASHING CHECK CONDUCTIVITY SENSOR**

- allows to detect impurities in the washing lye
- activates further rinse phases for achieving the perfect rinse



**Machines**

GW6010L	GW6010LX	GW6010ML	GW6010MLX	GW6010T
GW6010TM	GW7010	GW7010M	GW7010T	GW7010TM
GW7015	GW7015T			

**IC6290** (art. 904554)

**QUALITY WASHING CHECK CONDUCTIVITY SENSOR**

- allows to detect impurities in the washing lye
- activates further rinse phases for achieving the perfect rinse



**Machines**

GW6290DS

**IRCP** (art. 903000)

**CEILING SPRAY ARM**

- additional ceiling spray arm
- improves the washing efficiency of complex glassware shapes
- made of stainless steel AISI 316L
- not compatible with the trolleys CSK6B, CSK8B, LM40S, LB16S, CSK4BMIC



**Machines**

GW0260-00-00J	GW0260-S0-00J	GW1260-00-000	GW1260-S0-000
GW1260-SC-000	GW3060BX	GW3060BXA	GW3060BXC
GW3060BXCA	GW3060P	GW3060PA	GW3060S
GW3060SA	GW3060SC	GW3060SCA	GW4090

GW4090A	GW4090BXC	GW4090BXCA	GW4090BXP
GW4090C	GW4090CA	GW4090P	GW4090PA
GW4260-00-000	GW4260-S0-000	GW4260-S0-060	GW4260-SC-000
GW4290-00-000	GW4290-S0-000		

**KITSUPPR (art. 904886)**

**KIT WHEELS SUPPORT FOR THE UPPER BASKET IN UNDERCOUNTER**

- Kit of reinforced support for the wheels of the upper basket into the undercounter machines
- compatible with all the undercounter machine models provided of internal upper level wheels
- better and quicker replacement "plug and play" of internal wheels in phase of routine maintenance of the machine



**Machines**

GW0260-00-00J	GW0260-S0-00J	GW1260-00-000	GW1260-S0-000
GW1260-SC-000	GW3060BX	GW3060BXA	GW3060BXC
GW3060BXCA	GW3060P	GW3060PA	GW3060S
GW3060SA	GW3060SC	GW3060SCA	GW4090
GW4090A	GW4090BXC	GW4090BXCA	GW4090BXP
GW4090C	GW4090CA	GW4090P	GW4090PA
GW4260-00-000	GW4260-S0-000	GW4260-S0-060	GW4260-SC-000
GW4290-00-000	GW4290-S0-000		

**LED6010 (art. 903409)**

**LED LIGHTING FOR WASHING CHAMBER**

- n. 4 LED lights ensuring complete chamber illumination
- allows chamber inspection during the washing cycle
- lights can be replaced from inside the chamber



**Machines**

GW6010L	GW6010LX	GW6010ML	GW6010MLX	GW7010
GW7010M	GW7015			

## **STEAMCO15** (art. 904267)

---

### **ECOSLIM STEAM CONDENSER**

- high efficiency condenser
- cold water spray
- to be used for eliminating exhaust air vapour
- avoids connection with exhaust air duct

#### **Machines**

GW7015                      GW7015T

## **STEAMCO60** (art. 903466)

---

### **ECOSLIM STEAM CONDENSER**

- high efficiency condenser
- cold water spray
- to be used for eliminating exhaust air vapour
- avoids connection with exhaust air duct

#### **Machines**

GW6010L                      GW6010LX                      GW6010ML                      GW6010MLX                      GW6010T  
GW6010TM                      GW7010                      GW7010M                      GW7010T                      GW7010TM

## **STEAMFEED15** (art. 904313)

---

### **HYBRID STEAM HEATING FOR CHAMBER**

- steam heating element for washing chamber
- efficiently reduces electrical consumption
- hybrid system electric + steam
- provided with condensed steam "loop"
- made of stainless steel

#### **Machines**

GW7015                      GW7015T

## STEAMFEED15T (art. 904314)

### HYBRID STEAM HEATING FOR CHAMBER AND PRE-HEATING TANKS

- steam heating element for washing chamber and pre-heating tanks
- efficiently reduces electrical consumption
- hybrid system electric + steam
- provided with condensed steam "loop"
- made of stainless steel

#### Machines

GW7015T

## STEAMFEED70 (art. 904016)

### HYBRID STEAM HEATING

- steam heating element for washing chamber
- reduces efficiently electrical consumption
- hybrid system electric + steam
- provided with condensed steam "loop"
- made of stainless steel

#### Machines

GW7010

GW7010M

GW7010T

GW7010TM

## T4260AC (art. 905442)

### HOT WATER SUPPLY SYSTEM

- optional system for hot water machine supplying
- provided with water hose, Aquastop safety system included



#### Machines

GW0260-00-00J

GW0260-S0-00J

GW1260-00-000

GW1260-S0-000

GW1260-SC-000

GW4260-00-000

GW4260-S0-000

GW4260-S0-060

GW4260-SC-000

GW4290-00-000

GW4290-S0-000

**TOP45I** (art. 902032)

**BUILT-IN TOP 45 CM**

- for built-in mounting solution
- reduces machine height at 830 mm



**Machines**

GW2145                      GW2145A

**TOP60I** (art. 902033)

**BUILT-IN TOP 60 CM**

- for built-in mounting solution
- reduces machine height at 830 mm



**Machines**

GW0260-00-00J	GW0260-S0-00J	GW1260-00-000	GW1260-S0-000
GW1260-SC-000	GW3060BX	GW3060BXA	GW3060BXC
GW3060BXCA	GW3060P	GW3060PA	GW3060S
GW3060SA	GW3060SC	GW3060SCA	GW4260-00-000
GW4260-S0-000	GW4260-S0-060	GW4260-SC-000	GW4290-00-000
GW4290-S0-000			

**TOP90I** (art. 902034)

**BUILT-IN TOP 90 CM**

- for built-in mounting solution
- reduces machine height at 830 mm



**Machines**

GW4090	GW4090A	GW4090BXC	GW4090BXCA	GW4090BXP
GW4090C	GW4090CA	GW4090P	GW4090PA	

## TOP90IC (art. 903523)

### BUILT-IN TOP 90 CM

- for built-in mounting solution
- reduces machine height at 830 mm



### Machines

GW4290-00-000

GW4290-SC-000

## WD-CHECK15 (art. 905102)

### CIRCUIT FOR WATER SAMPLE WITHDRAWAL

- Special circuit for sampling, for taking water samples from the wash tank
- The sampling tap is positioned on the loading side, accessible by opening the door of the tank compartment
- Available for 7015 series machines



### Machines

GW7015

GW7015T

## WD-CODESCAN (art. 905103)

### BARCODE READER

- CODE128 barcode reader
- connected to the machine by cable
- positioned on the loading side
- available for 7010-7015 series machines - not compatible with loading side printer



### Machines

GW7010

GW7010M

GW7010T

GW7010TM

GW7015

GW7015T

**WD-DRAIN15** (art. 905101)

**WASHING PUMP DRAIN CIRCUIT**

- Draining circuit of the washing pumps
- It acts at the end of each phase of the washing cycle
- Guarantees complete emptying of the pump
- Available for 7015 series machines



**Machines**

GW7015                      GW7015T

**WD-EC-1** (art. 905558)

**SINGLE-PHASE ELECTRIC CONNECTION**

- for converting the machine from the standard configuration (3/N/PE 400 V – 50 or 60 Hz) to the single-phase configuration 1/N/PE 230 V – 50 or 60 Hz
- NOTE: an authorized engineer is required for installation outside the factory

**Machines**

GW0260-00-00J      GW0260-S0-00J      GW1260-00-000      GW1260-S0-000      GW1260-SC-000  
 GW4260-00-000      GW4260-S0-000      GW4260-S0-060      GW4260-SC-000

**WD-EC-3** (art. 905559)

**THREE-PHASE WITHOUT NEUTRAL ELECTRIC CONNECTION**

- for converting the machine from the standard configuration (3/N/PE 400 V – 50 or 60 Hz) to the three-phase without neutral 3/PE 230 V - 50 or 60 Hz
- NOTE: an authorized engineer is required for installation outside the factory

**Machines**

GW0260-00-00J      GW0260-S0-00J      GW1260-00-000      GW1260-S0-000      GW1260-SC-000  
 GW4260-00-000      GW4260-S0-000      GW4260-S0-060      GW4260-SC-000

## WD-EC-1UK (art. 905726)

---

### SINGLE-PHASE ELECTRIC CONNECTION WITH UK PLUG

- for converting the machine from the standard configuration (3/N/PE 400 V – 50 Hz) to the single-phase configuration, with UK plug and 13A protection
- NOTE: installation outside the factory requires the intervention of an authorized technician

#### Machines

GW4260-00-000      GW4260-0C-000      GW4260-S0-000      GW4260-S0-060      GW4260-SC-000  
GW4260-SC-00G

## WD-EC4290-1 (art. 905614)

---

### SINGLE-PHASE ELECTRIC CONNECTION

- Kit of components needed to convert the machine from the standard configuration (3/N/PE 400 V – 50 or 60 Hz)
- to single-phase configuration 1/N/PE 230 V – 50 or 60 Hz
- the components are SPECIFIC for GW and WD4290
- NOTE: A technician is required for installation

#### Machines

GW4290-00-000      GW4290-S0-000      GW4290-SC-000

## WD-EC4290-3 (art. 905615)

---

### THREE-PHASE WITHOUT NEUTRAL ELECTRIC CONNECTION

- kit of components needed to convert the machine from the standard configuration (3/N/PE 400 V – 50 or 60 Hz)
- to three-phase configuration without neutral 3/PE 230 V - 50 or 60 Hz
- the components are SPECIFIC for GW and WD4290
- NOTE: A technician is required for installation


#### Machines

GW4290-00-000      GW4290-S0-000      GW4290-SC-000

## WD-EMS (art. 905104)

### ABSORBED POWER MANAGEMENT SYSTEM

- Application EMS - Energy Management System - for managing the power absorbed by machines installed in battery in the sterilization unit or big laboratories
- It displays the progress of the cycles in execution of the connected machines
- Available for 6010/7010/7015 series machines

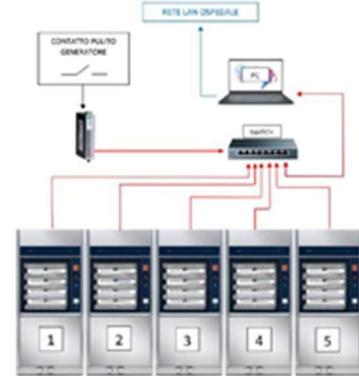

### Machines

GW6010L	GW6010LX	GW6010ML	GW6010MLX	GW6010T
GW6010TM	GW7010	GW7010M	GW7010T	GW7010TM
GW7015	GW7015T			

## WD-LAN4290 (art. 905443)

### COMMUNICATION LAN CARD

- data communication card
- for machines of the 4290 product family
- equipped with RS232 serial port for connection to an external printer
- equipped with a LAN port for ethernet connection which allows communication and monitoring of the machine from any access point in the local network
- to be set using the instructions in the installation manual



## WD-SIGNBOX15 (art. 905135)

### FREE CONTACTS BOX FOR REMOTE SIGNAL

- Box containing 3 free contacts
- 1 general alarm contact
- 1 warning contact
- 1 contact that identifies if the cycle is in progress
- Positioned on the top of the machine
- Available for 7010 and 7015 families



### Machines

GW7010	GW7010M	GW7010T	GW7010TM	GW7015
GW7015T				

## WD-SIGNBOXIN (art. 905184)

### FREE CONTACTS BOX FROM EXTERNAL INPUT

- Box with 1 external clean contact supplied by electric generator, for 7015 and 7010 families
- If contact open -> generator off
- If contact closed -> generator on
- To be connected to WD-EMS to command its automatic start



### Machines

GW7010	GW7010M	GW7010T	GW7010TM	GW7015
GW7015T				

## WD-VDS (art. 902592)

### DUAL DRAIN VALVE

- separates and channels polluted drain water used in the first washing cycle for final rinsing water



### Machines

GW3060BX	GW3060BXA	GW3060BXC	GW3060BXCA	GW3060P
GW3060PA	GW3060S	GW3060SA	GW3060SC	GW3060SCA
GW4090	GW4090A	GW4090BXC	GW4090BXCA	GW4090BXP
GW4090C	GW4090CA	GW4090P	GW4090PA	GW6010L
GW6010LX	GW6290DS			

## WD-VDS6010 (art. 903671)

### DUAL DRAIN VALVE

- separates and channels polluted drain water used in the first washing cycle for final rinsing water
- to be necessarily combined with wall drain pump (art. DRAINP6010)



### Machines

GW6010ML	GW6010MLX	GW6010T	GW6010TM	GW7010
GW7010M	GW7010T	GW7010TM	GW7015	GW7015T

## WD-WS6010 (art. 905553)

---

### 6010 WATER SOFTENER

- onboard water softener, with automatic regeneration managed by the processor
- salt lack sensor with visual signal on the display
- wires and brackets customized for GW and WD6010 family
- optional installed at factory

#### Machines

GW6010L	GW6010LX	GW6010ML	GW6010MLX	GW6010T
GW6010TM				

## WD-WS7010 (art. 905553)

---

### 7010 WATER SOFTENER

- onboard water softener, with automatic regeneration managed by the processor
- salt lack sensor with visual signal on the display
- wires and brackets customized for GW and WD7010 family
- optional installed at factory

#### Machines

GW7010	GW7010M	GW7010T	GW7010TM
--------	---------	---------	----------

## WD-IP216010 (art. 905707)

---

### IP21 PROTECTION LEVEL

- allows the achievement of IP21 protection level: protection against dripping water
- for WD6010 version
- factory optional



#### Machines

GW6010L	GW6010LX	GW6010ML	GW6010MLX
---------	----------	----------	-----------

## **WD-IP216010T** (art. 905708)

---

### **IP21 PROTECTION LEVEL**

- allows the achievement of IP21 protection level: protection against dripping water
- for WD6010T version
- factory optional



### **Machines**

GW6010T                      GW6010TM

## **WD-IP217010** (art. 905709)

---

### **IP21 PROTECTION LEVEL**

- allows the achievement of IP21 protection level: protection against dripping water
- for 7010/7010T versions
- factory optional



### **Machines**

GW7010                      GW7010M

## Plinths and support frames

### B6040L (art. 902487)

#### STAINLESS STEEL PLINTH FOR 60 CM MACHINES WITH LOCK

- dimensions LxDxH: 600 x 530 x 400 mm
- increases the machine working level to 70 cm from the floor
- detergent cabinet with front door
- equipped with adjustable feet
- detergent cabinet with anti-drop tray
- front door with lock
- made of stainless steel



#### Machines

GW0260-00-00J	GW0260-S0-00J	GW1260-00-000	GW1260-S0-000	GW1260-SC-000
GW3060BX	GW3060BXA	GW3060BXC	GW3060BXCA	GW3060P
GW3060PA	GW3060S	GW3060SA	GW3060SC	GW3060SCA
GW4260-00-000	GW4260-S0-000	GW4260-S0-060	GW4260-SC-000	GW4290-00-000
GW4290-S0-000				

### B9040L (art. 902813)

#### STAINLESS STEEL PLINTH FOR 90 CM MACHINES WITH LOCK

- dimensions LxDxH: 900 x 530 x 400 mm
- increases the machine working level to 70 cm from the floor
- detergent cabinet with front door
- equipped with adjustable feet
- detergent cabinet with anti-drop tray
- front door with lock
- made of stainless steel



#### Machines

GW4090	GW4090A	GW4090BXC	GW4090BXCA	GW4090BXP
GW4090C	GW4090CA	GW4090P	GW4090PA	GW4290-00-000
GW4290-S0-000	GW4290-SC-000			